Product data sheet



Hinged support rail Mono WARM TOUCH 950.50.11050







Hinged support rail Mono WARM TOUCH 950.50.11050


- with a horizontal, round gripping level with diameter 33 mm
- · for holding on to and providing support
- 600 mm projection, wall mounting plate 181 mm high and 130 mm wide
- · can be folded upwards and, folded downwards with braking mechanism
- made of high-quality chrome-look coated polyamide that is warm to touch
- with continuous anti-corrosion steel core and integrated wall mounting plate made of steel with cover made of chrome-look coated polyamide
- · stable three-point fixing
- · for wall mounting with concealed, non-corrosive and tested fixing materials for different wall systems
- with toilet roll holder 950.50.01050, toilet roll holder with WC flushing mechanism and function button (radiocontrolled) 950.50.06290 (right) and 950.50.06390 (left) for retrofitting
- CE marking: Medical devices Class 1 according to Directive 93/42 / EEC
- · fulfils the DIN 18040 requirements







Dimensions in mm

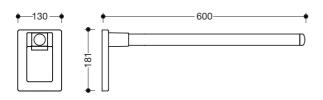

Surface

WARM TOUCH

Chrome-look coated polyamide

Technical information subject to alteration, 15.06.2015