



M9S5E




Hinged support rail Duo WARM TOUCH 950.50.12050

- with two gripping levels arranged above each other
 bottom rail is merged with the top rail at an angle of 135°
- top rail (diameter 33 mm) for holding on to and providing support
- bottom rail (diameter 25 mm) starts at a distance of 250 mm from the front edge of the top rail and is used for holding on to and for turning in from the side
- 700 mm projection, wall mounting plate 181 mm high and 130 mm wide
- \cdot can be folded upwards and, folded downwards with braking mechanism
- \cdot made of high-quality chrome-look coated polyamide that is warm to touch
- with continuous anti-corrosion steel core and integrated wall mounting plate made of steel with cover made of chrome-look coated polyamide
- · stable three-point fixing
- · for wall mounting with concealed, non-corrosive and tested fixing materials for different wall systems
- with toilet roll holder 950.50.01050, toilet roll holder with WC flushing mechanism and function button (radiocontrolled) 950.50.06290 (right), 950.50.06390 (left), and floor support 950.50.02092 for retrofitting
- CE marking: Medical devices Class 1 according to Directive 93/42 / EEC
- · fulfils the ÖNORM B1600/1601 and SIA 500 requirements



Dimensions in mm

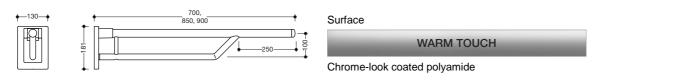

Technical information subject to alteration, 15.06.2015

HEWI Heinrich Wilke GmbH Postfach 1260 D-34442 Bad Arolsen Telefon: +49 5691 82-366 Telefax: +49 5691 82-99366 international@hewi.com www.hewi.com

HEWI (UK) Limited Holm Oak Barn, Beluncle Halt Stoke Road Hoo Kent ME3 9NT Phone +44 1634 258 200 Fax +44 1634 250 099 info@hewi.co.uk www.hewi.co.uk